

KB03007 Thiol Quantification Assay Kit

96 well plate 100/200/400 tests





Table of contents

١.	General information	l
2.	Technical specifications	2
3.	Materials and storage	3
4.	Introduction	2
5.	Assay principle	2
6.	Assay preparation	5
7.	Sample preparation	7
8.	Assay protocol	8
9.	Data analysis	10
10.	Troubleshooting	12
11.	Additional information	14
12.	Related products	14
13.	Warranties and limitation of liability	1.5



1. General information

PRECAUTIONS

Please read this manual carefully before beginning the assay.

This product is designed for **research use only**. It is not approved for human or animal use or clinical diagnosis. All chemicals should be handled with care and in accordance with laboratory safety practices. It is recommended to use basic Personal Protective Equipment.

Do not use after the expiration date stated on the packaging.

Do not mix or substitute reagents or materials from other kit batches or vendors.

For the **material safety data sheet** (MSDS) please contact us at info@bioquochem.com

TECHNICAL RECOMMENDATIONS

Store reagents as indicated in Materials and storage section.

Be sure to keep the bottle capped when not in use.

Let the components reach room temperature (RT) before use.

Immediately before use, gently invert and rotate reagent bottles several times to mix the contents thoroughly.

Avoid foaming or bubbles when mixing or reconstituting components.

Avoid cross contamination of samples or reagents by changing pipette tips between sample, standard and reagent additions.

Be sure to use the optimal microplate for the assay. Flat bottom transparent microplates for UV/VIS applications, and black microplates for fluorescence measurements.



2. Technical specifications

Available sizes

100/200/400 tests

Required sample volume

20 µL/test

Compatible samples

Biological fluids, cell lysates and tissue homogenates

Type of detection

Colorimetric (412 nm)



3. Materials and storage

MATERIALS SUPPLIED

Item	No. Tests	Units	Storage
	100	1	
Reagent A	200	1	4 °C
	400	2	
	100	1	
Reagent B	200	1	4 °C
	400	2	
	100	1	
Reagent C	200	2	4 °C
	400	4	
	100	1	
Reagent D	200	2	4 °C
	400	4	
	100	1	
Reagent E	200	1	4 °C
	400	1	
	100	1	
Standard	200	2	4 °C
	400	4	
	100	1	
Transparent 96-Well Microplate	200	2	RT
	400	4	

MATERIALS NEEDED BUT NOT SUPPLIED

- o Double distilled water (ddH2O) as Milli-Q Ultrapure Water
- o Labware materials (micropipettes, tubes, stirring/mixing equipment)
- Colorimetric microplate reader equipped with filter for OD 412 nm

STORAGE CONDITIONS

On receipt, store kit components as indicated above. Under these conditions, the reagents are stable in the original packaging until the expiration date indicated on the outside of the box. After reconstitution, standard solutions are unstable in the presence of oxygen. Prepare a fresh set of standards for every use.



4. Introduction

Protein (cysteine residues) and non-protein (e.g. homocysteine and glutathione) thiol groups (R-SH) are central components of the extracellular antioxidant defense. High concentrations of free thiols are indicative of a favorable redox status, while low levels of free thiols are associated with different oxidative stress related diseases. Reduced levels of free thiols arise from its rapid oxidation to disulfides in the presence of high amounts of reactive oxygen species (ROS).

BQC Thiol Quantification Assay Kit is a selective, easy, and quick assay for measuring free thiols and disulfides groups in biological samples.

5. Assay principle

BQC Thiol Quantification Kit is based on the Ellman's method. The reaction between the Ellman's Reagent (DTNB: 5,5'-dithiobis (2-nitrobenzoic acid)) and a thiol group (R-SH) leads to the formation of 2-nitro-5- thiobenzoate (TNB-) that ionizes to the yellow colored (412 nm) dianion TNB²⁻ at alkaline pH. The intensity of the color produced is proportional to thiol concentration.

R-SH + HO
$$O_2N$$
 OH R-S-S-R + O_2N NO O_2 R-S-S-R + O_2N Thiol groups DTNB Disulfide TNB2-Yellow Color O_2N O_2

Reaction of DTNB with thiol group

In addition to free thiol determination, this kit also allows for the quantification of total thiols by reducing the disulfides groups present in a sample with NaBH₄.

Schematic determination of free and total thiols



6. Assay preparation

REAGENT PREPARATION

All assay reagents not listed below are ready to use as supplied. Allow the reagents to reach room temperature before use.

- **R.B. Working Solution (only for total thiol determination)**: Add 2 mL of Reagent E and 2 mL of ddH₂O to each vial of Reagent B and mix well.
 - CAUTION: R.B. Working Solution must be freshly prepared and discarded after use
- **R.D. Working Solution**: Add 2.5 mL of Reagent E to each vial of Reagent D and mix well.
 - CAUTION: R.D. Working Solution must be freshly prepared and discarded after use

Standard Solution (Glutathione, GSH): Add 2 mL of ddH₂O to the Standard vial and mix well. Dilute this standard solution 1:10 with ddH₂O (e.g. 100 μ L of standard solution with 900 μ L of ddH₂O). Use this 1:10 diluted Standard solution to prepare the standard curve.

STANDARD CALIBRATION

Prepare Glutathione (GSH) standards for the calibration curve from the diluted Standard solution according to the following Table. Prepare the standards immediately prior to each assay. Vortex tubes thoroughly. Discard standard solutions after use.

Standard	Standard solution 1:10 diluted (µL)	ddH ₂ O (μL)	GSH (mM)
Std 1 (Reagent Blank)	0	100	0
Std 2	20	80	0.2
Std 3	40	60	0.4
Std 4	60	40	0.6
Std 5	80	20	0.8
Std 6	100	0	1.0



PLATE SET UP

BQC recommends running the standards and samples at least in duplicate (triplicate recommended). There is no specific pattern for using the wells on the plate. A proposed layout of standards (Std) and samples (S) to be measured in duplicate is shown below.

NOTE: If sample blanks are included in the assay, it is necessary to reserve some wells of the plate for these blanks

Q	1	2	3	4	5	6	7	8	9	10	11	12
Α	Std 1	Std 1	S3	S3	S11	S11	S19	S19	S27	S27	\$35	\$35
В	Std 2	Std 2	S4	S4	\$12	S12	S20	S20	S28	S28	\$36	\$36
С	Std 3	Std 3	\$5	S5	\$13	\$13	S21	S21	S29	S29	\$37	S37
D	Std 4	Std 4	S6	S6	\$14	\$14	S22	S22	S30	\$30	\$38	\$38
E	Std 5	Std 5	S7	S7	\$15	\$15	S23	S23	S31	S31	S39	S39
F	Std 6	Std 6	S8	S8	\$16	\$16	S24	S24	S32	S32	S40	\$40
G	S1	S 1	S9	S9	\$17	S17	S25	\$25	\$33	\$33	S41	S41
Н	S2	S2	S10	\$10	\$18	\$18	S26	S26	\$34	\$34	S42	S42

Example of plate layout for the Thiol Assay Kit



7. Sample preparation

The following sample preparation protocols are intended as a guide only. The optimal conditions for sample preparation must be determined by the end user. It is recommended to use fresh samples. If it is not possible, aliquot and store samples appropriately with minimal freeze/thawing.

Thiol Assay Kit can be used to determine thiols and disulfides in a wide variety of biological samples like serum, plasma, urine, cell lysate or tissue homogenate.

Biological samples. Biological samples like plasma, serum or urine can be directly measured with appropriate dilutions.

Cell lysate. Wash cells with ice-cold buffer before lysis. Lyse cells by sonication or freeze-thaw cycles. Centrifuge cell lysis suspension at 10000 x g for 15 minutes at 4 $^{\circ}$ C and collect the supernatant. It is recommended to use lysates from $2 \cdot 10^6$ cells.

Tissue Homogenates. Dissect the tissue of interest and place it on a homogenizer tube with an appropriate amount of an ice-cold buffer (e.g. 1 g tissue per 1-2 mL buffer). Homogenize the tissue and then centrifuge the homogenate at 10000 x g for 15 minutes at 4 °C. Collect the supernatant.

Reagents and materials required for sample preparation are not supplied with the kit. Before doing sample preparation, consider the volume of sample required per test; see **Technical specifications** section.

Make sure that interfering substances present in the sample do not give a significant background. Run proper blanks as necessary (e.g. sample blank should be always evaluated when working with highly colored samples). It is recommended to assay different sample dilutions to ensure the values fall within the standard curve.



8. Assay protocol

Prepare and mix all reagents thoroughly before use. Each standard, sample or blank should be assayed at least in duplicate.

It is possible to calculate both **free thiols** (naturally reduced) and **total thiols** (naturally + chemically reduced) in each sample.

A. Free thiols determination:

1	Set up the plate design
2	Add 20 µL of standard or sample in each well
3	Add 20 μL of Reagent A in each well and wait 10 minutes
4	Add 220 µL of Reagent C in each well
5	Add 20 μL of R.D. Working Solution in each well and wait 10 minutes
6	Read the absorbance of all wells at 412 nm in end point mode at RT



B. Total thiols determination

1	Set up the plate design
2	Add 20 µL of standard or sample in each well
3	Add 20 µL of R.B. Working Solution in each well and wait 10 minutes
4	Add 220 µL of Reagent C in each well
5	Add 20 μL of R.D. Working Solution in each well and wait 10 minutes
6	Read the absorbance of all wells at 412 nm in end point mode at RT


If you need to **adapt this kit** for another form of the assay (for example cuvette), **contact us at** info@bioquochem.com



9. Data analysis

ANALYSIS OF THE STANDARDS

- Calculate the average absorbance of all the standards.
- Subtract the average absorbance of the reagent blank (Std 1) from the average absorbance of all the standards to obtain the blankcorrected absorbance of the standards.
- Create a standard curve by plotting the blank-corrected absorbance of the standards as a function of the standard concentration (see STANDARD CALIBRATION section). A typical standard curve (y=slope·x ± intercept) for this assay is shown below.

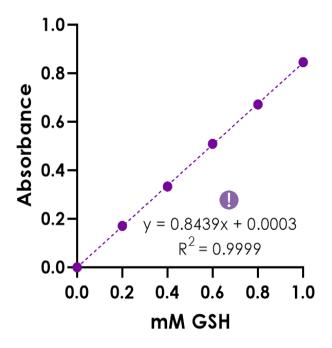

GSH standard curve with Thiol Assay Kit

This standard curve is an example of the data typically obtained with this kit. DO NOT USE this standard curve to calculate the thiol content of your samples. A new standard curve must be performed by the end user.



ANALYSIS OF THE SAMPLES

Both for free and total thiols

- Calculate the average absorbance of the samples.
- Subtract the average absorbance of the reagent blank (Std 1) from the average absorbance of each sample to obtain the blank-corrected absorbance of the samples (A_s).
- Calculate the GSH concentration of the samples using the following equation. Slope and intercept values are obtained from the standard curve.

GSH (mM) =
$$\left(\frac{A_S - intercept}{slope}\right)$$

When working with diluted samples the GSH concentration values obtained must be multiplied by the dilution factor to obtain the thiol content of the undiluted sample.

The results obtained depend on the assay format:

Free thiols (naturally reduced): concentration value obtained for the **free thiols assay**.

Total thiols (chemically + naturally reduced): concentration value obtained for the **total thiols assay**.

Disulfides: concentration obtained by **subtracting** the concentration calculated for **free thiols from the total thiols** concentration.



10. Troubleshooting

This troubleshooting table provides potential sources and solutions for common problems observed with BQC Assay Kits. **The problems listed below could occur when using any BQC Assay Kit**. They are not specific for this assay kit.

Problem	Possible Cause	Recommended Solution
Wells have color but there is no reading	Plate read at incorrect wavelength	Check the wavelength used in the assay
	Incorrect microplate	Use the correct microplate for your application UV/Vis: transparent Fluorescence: black wells/transparent bottom
	Pipetting errors in preparation of standards	Avoid pipetting small volumes (<5 µL) Be careful not to splash from well to well
	Air bubbles formed in well(s)	Use reverse pipetting technique
Standard readings do not	Standard stock is at incorrect concentration	Always refer to dilutions described in Assay preparation
follow a linear pattern	Improperly thawed reagents	Thaw all components completely and mix well before use
	Use of improperly stored reagents	Store the components appropriately Use fresh components from the standard curve
	Incorrect incubation times or temperatures	Refer to Assay protocol
Dispersion of standard and sample readings	Pipetting errors	Avoid pipetting small volumes (<5 µL) Be careful not to splash from well to well
	Air bubbles formed in well(s)	Use reverse pipetting technique



Problem	Possible Cause	Recommended Solution
	Samples contain interfering substances	Dilute sample further (if possible)
Sample erratic	Inappropriately stored samples or samples used after multiple freeze-thaw cycles	Use fresh samples or store appropriately until use
values	Samples not deproteinized	Use an appropriate deproteinization protocol
	Cells/Tissue samples not homogenized completely	Repeat the sample homogenization
	Inappropriate sample dilution buffer	Refer to Assay preparation
Sample reading fall outside the detection range	Samples are too diluted/concentrated No analyte/activity is observed in the sample	Re-assay using different sample dilutions

STILL HAVING PROBLEMS?

Please, contact BQC if you have any further questions, our team will be pleased to help you:

Phone	+ 34 985 26 92 92
E-mail	info@bioquochem.com
Business hours	Monday-Thursday: 8.30 to 17.00 (CEST) Friday: 8.00 to 15.00 (CEST)



11. Additional information

Thiol Quantification Assay Kit is a quick (< 30 minutes) and highly reproducible assay (RSD < 5 %) for free thiols and disulfide groups quantification in biological samples.

Cyanide, hydrogen sulphite and sulphide and some saturated and unsaturated aliphatic aldehydes have been reported to interfere with this assay.

If unexpected results are obtained running your samples, please contact us at info@bioquochem.com

12. Related products

More products available on bioquochem.com

Reference	Product
KB03016	MDA-TBARs Assay Kit
KB03032	Xanthine Oxidase Activity Assay Kit
KF01004	ORAC Antioxidant Capacity Assay



13. Warranties and limitation of liability

BQC shall not in any event be liable for incidental, consequential or special damages of any kind resulting from any use or failure of the products, even if BQC has been advised of the possibility of such damage including, without limitation, liability for loss of use, loss of work in progress, downtime, loss of revenue or profits, failure to realize savings, loss of products of buyer or other use or any liability of buyer to a third party on account of such loss, or for any labor or any other expense, damage or loss occasioned by such product including personal injury or property damage is caused by BQC's gross negligence. Any and all liability of BQC hereunder shall be limited to the amounts paid by the buyer for the product.

Buyer's exclusive remedy and BQC's sole liability hereunder shall be limited to a refund of the purchase price, or the replacement of all material that does not meet our specifications.

Said refund or replacement is conditioned on buyer giving written notice to BQC within 30 days of shipment.

Expiration date: 1 year from the date of fabrication. Expiration date is indicated on the outside of the box.

For further details, please refer to our website bioquochem.com



Edificio CEEI | Parque Tecnológico de Asturias,

33428 Llanera, Asturias

Info@bioquochem.com



www.bioquochem.com