

Colony Counter



Features

4 panel displays

4 panel display offers convenience of showing the present and previous three count values.

Multi-counting average

Galaxy 330 can make average count from multiple counting.

Exclusive BACK button

The BACK button offers the convenience of deducting data when over counting.

Suitable for Petri-dishes in various diameters

Galaxy 330's adjustable dish holder can accept most Petri-dishes in a wide range of diameters from 90 to 150 mm or 60 mm.

Pressure sensor design

Pressure sensor design can allow user to operate with any marking pens. No worry about sample spoiled or contaminated.

International certification

CE certification

Two-Year Warranty

Applications

Microbiology colony counting

Total bacteria counting

Total coliform counting

Ordering Information

175330-01(02)

Galaxy 330, Colony Counter with AC100-240V adaptor, US plug (EU plug)


175331-01(02)

Galaxy 330, Colony Counter with software, AC100-240V adaptor, US plug (EU plug)

175330-12

Software Package for Galaxy330 (Include a 8G USB flash drive and connecting cable)

Galaxy 330

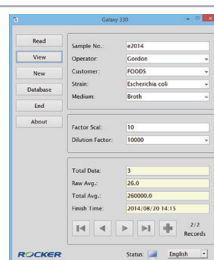

Galaxy 330, PC Software (optional)

In addition to recording and tracing data, the Galaxy 330 PC software can help user add, delete, search, print out and save data as Excel format (csv file).

Specification

Model	Galaxy 330
Display	3-digit (0-999) LED display, 4 panels
Magnification	2~3 X
Petri dish diameter	60 mm or 90-150 mm (adjustable)
Memory	100
Average	Yes
Zero reset	Yes
Count back	Yes
Illumination	LED (adjustable)
Pressure sensitivity	Adjustable
Acoustical counting signal	Adjustable
Output	USB
Dimension (LxWxH)	35.1 x 30.2 x 37 cm
Net weight	5.6 Kg