

# Stainless Steel Spin-Lock Design Multi-Branch Manifold



VIDEO

2  
YEAR  
WARRANTY

## Features

### Quick-assembly and disassembly

Each component can be easily removed for cleaning and autoclaving. All wetted parts of MultiVac 300 / 600 manifold are made of SS316 stainless steel, which has a better chemical resistance.

### Convertible funnel adaptor

Convertible funnel adaptor design allows user to fit the desired funnel (reusable or disposable) easily.

### Innovative vent valve design

An adjustable vent valve integrated with the pipe can eliminate water hammer phenomenon and help drain out residual waste water.

### Quick-Fit hose barb

Detachable hose barb can be connected to either side of the manifold quickly towards the vacuum source.

### Individual control valve

Each branch permits individual control.

### Compatible with different sizes of funnels

MultiVac 300 - MB / 600 - MB manifolds can accept Rocker standard SS 100, 300, 500 ml funnel at option.

### Two-Year Warranty

## Applications

- Microbiological test
- Chemical experiment
- Vacuum filtration

## Ordering Information

### 180300-02

MultiVac 300 - MB, 3 Places SS Manifold with Rocker SS funnel adaptor

### 180600-02


MultiVac 600 - MB, 6 Places SS Manifold with Rocker SS funnel adaptor

**167103-61** 100 ml SS Funnel

**167103-63** 300 ml SS Funnel

**167103-65** 500 ml SS Funnel

### MultiVac 300 - MB

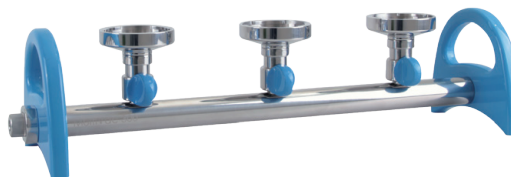


Adjustable vent valve



Quick-Fit hose barb

### MultiVac 600 - MB

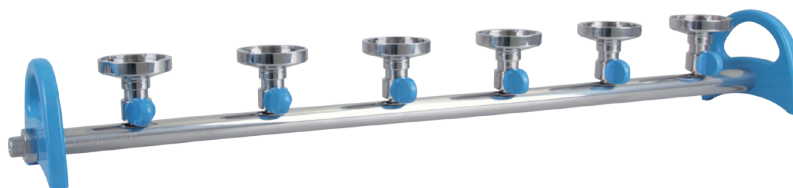

## Specification

Model	MultiVac 300 - MB	MultiVac 600 - MB
Number of branch	3	6
Filter diameter	47 mm	47 mm
Funnel capacity	100 ml / 300 ml / 500 ml	100 ml / 300 ml / 500 ml
Effective filtration area	9.6 cm <sup>2</sup>	9.6 cm <sup>2</sup>
Hose barb	5/16 inch (8 mm)	5/16 inch (8 mm)
Dimension (LxWxH)	48 x 13 x 11 cm	85 x 13 x 11 cm
Part	Material	Material
Body	SS316	SS316
Control valve	SS316	SS316
Adaptor	SS316	SS316
Membrane support	SS316	SS316
Handle / Knob	Aluminum	Aluminum
Hose barb / Vent valve	SS316	SS316

# Stainless Steel Spin-Lock Design Multi-Branch Manifold



## Features

### Quick-assembly and disassembly

Each component can be easily removed for cleaning and autoclaving. All wetted parts of MultiVac 300 / 600 manifold are made of SS316 stainless steel, which has a better chemical resistance.

### Quick-Fit hose barb

Detachable hose barb can be connected to either side of the manifold quickly towards the vacuum source.

### Innovative vent valve design

An adjustable vent valve integrated with the pipe can eliminate water hammer phenomenon and help drain out residual waste water.

### Individual control valve

Each branch permits individual control.

### Compatible with Rocker standard 100 ml SS funnels

### Two-Year Warranty

## Applications

- Microbiological test
- Chemical experiment
- Vacuum filtration

## Ordering Information

### 180301-02

MultiVac 301 - MB, 3 Places SS Manifold with Rocker SS funnel adaptor, 3 of 100 ml SS funnel

### 180601-02

MultiVac 601 - MB, 6 Places SS Manifold with Rocker SS funnel adaptor, 6 of 100 ml SS funnel

### MultiVac 301 - MB



Adjustable vent valve



Quick-Fit hose barb

### MultiVac 601 - MB



## Specification


Model	MultiVac 301 - MB	MultiVac 601 - MB
Number of branch	3	6
Filter diameter	47 mm / 50 mm	47 mm / 50 mm
Funnel capacity	100 ml	100 ml
Effective filtration area	9.6 cm <sup>2</sup>	9.6 cm <sup>2</sup>
Hose barb	5/16 inch (8 mm)	5/16 inch (8 mm)
Dimension (LxWxH)	48 x 13 x 11 cm	85 x 13 x 11 cm
Part	Material	Material
Body	SS316	SS316
Control valve	SS316	SS316
Adaptor	SS316	SS316
Membrane support	SS316	SS316
Handle / Knob	Aluminum	Aluminum
Hose barb / Vent valve	SS316	SS316

# How to order?

You can select Combination A, B or C to

**A**

Step 1. Choose the manifold

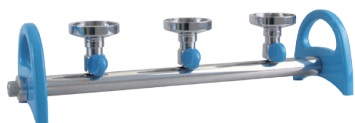


MultiVac 300 - MB  
3 Places Stainless Steel Manifold  
(180300-02)



MultiVac 600 - MB  
6 Places Stainless Steel Manifold  
(180600-02)

**B**

Step 1. Choose the manifold

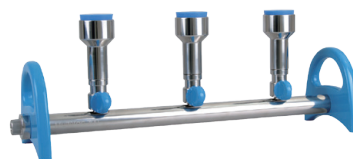

MultiVac 300 - MS  
3 Places Stainless Steel Manifold  
(180300-01)



MultiVac 600 - MS  
6 Places Stainless Steel Manifold  
(180600-01)

Step 2. Choose the SS funnel



Rocker, 100 ml  
SS Stainless Steel Funnel  
(167103-61)



Rocker, 300 ml  
SS Stainless Steel Funnel  
(167103-63)



Rocker, 500 ml  
SS Stainless Steel Funnel  
(167103-65)



Stainless Steel Cover  
for SS 100 ml funnel  
(167103-11)



Stainless Steel Cover  
for SS 300 ml funnel  
(167103-13)



Stainless Steel Cover  
for SS 500 ml funnel  
(167103-15)



Rocker, 300 ml  
MF 3, Magnetic Filter Holder  
(200300-01)



Rocker, 500 ml  
MF 5, Magnetic Filter Holder  
(200500-01)



Rocker, 300 ml  
VF6, Glass Filter Holder  
(167120-31)


# complete your order

C

Step 1. Choose the manifold



MultiVac 300  
3 Places Stainless Steel Manifold Base  
(180300-00)

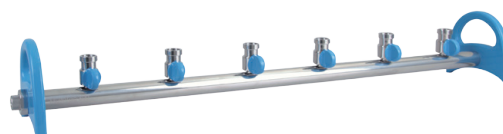

MultiVac 600  
6 Places Stainless Steel Manifold Base  
(180600-00)



Step 2. Choose the disposable funnel and adaptor



Sartorius, 100 ml  
disposable plastic funnel  
(16A07)



Millipore, 100 ml / 250 ml  
Microfil, disposable plastic funnel  
(MIACFD101 / MIACFD201)



Pall, 100 ml / 250 ml  
disposable plastic funnel  
(4870 / 4871)

+

+

+



SS316 Adaptor for Sartorius  
disposable plastic funnel  
(190100-74)



SS316 Adaptor for Millipore  
disposable plastic funnel  
(190100-73)



SS316 Adaptor for Pall  
disposable plastic funnel  
(190100-75)

\* SS316 adaptor is only available for above-mentioned disposable funnels.

## Customized Design

We can offer customized  
adaptor for your funnel

